Diamond Shield IE

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Mikroströme - EMS

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Firmware Version 2.1

Index

Introduction.....	29
Important Notice.....	31
Feature Identification.....	33
Bringing into service.....	36
Operation.....	36
Technical Data.....	42
For disposal.....	43
Declaration of Conformity.....	45

Introduction

The Diamond Shield IE Zapper

The Wellness Device with Microcurrents

- wide frequency range of 0.5 Hz to 2,5 MHz.
- small voltages (0.1 V) to high voltages (16 V)
- Impulse discharge
- Wobble functionality
- Modulation of the output frequencies
- Constant current control
- Running frequencies from special chip cards (optional accessory)

The operation can be done via a 9V battery or a rechargeable battery.

The device can be used to enhance the well-being.

An earthing cable and an enclosed power socket tester are available as optional accessories.

The good effect of the earthing (so called negative ions) occurs only if the power socket is properly earthed. Be sure to check this with the power socket tester.

Warning

The use of the earthing plug (or any other electrical appliances) to an incorrect power socket can be perilous.

Do not use the earthing cable during a thunderstorm.

Important Notice

Safety Instruction

- Use only the supplied grounding cable (optional accessory). Under no circumstances use a different cable for grounding
- Under any circumstances never insert, a plug into an electrical outlet
- Always insert the grounding cable to the Diamond Shield and then connect it to the power outlet.
- Never use this device to anyone who has a pacemaker.
- Never use this device in the eye area.
- Keep the Diamond Shield always clean and dry – wiping it with a damp cloth is permitted.

***Using of the device according to the operating instructions is essential.
DermaVit assumes no liability for misuse.***

Note

In case of feeling a tingle or burning on your skin during the usage please lower the intensity. The efficiency is independent from the above. When hypersensitive to metal, we gladly exchange the metal bracelets with cloth bracelets.

Unpacking

Before using this device, please check of the completeness of the supplied accessories:

- Connection cable for (hand-) electrodes or bracelets
- Grounding cable (optional accessory)
- 2 Bracelets
- 9V Battery

Storage and transportation

The storage of the equipment shall be always stored in dry, closed rooms. If the device is transported or stored in extreme temperatures, a period of at least 2 hours should be observed for acclimatization before handling.

Intended Usage

The device is used to increase the subjective well-being of people who are open to such an approach.

The devices are intended for use in clean, dry rooms. They cannot be used in extreme dust, damp, dusty air, danger of explosion or chemical vapors.

Protect the Diamond Shield from direct sunlight.

The permissible operating temperature ranges from +15°C to +30°C. During storage or transportation, the ambient temperature can be between -20°C and +70°C. If condensation has been formed during transportation or storage, the device must acclimate for about 2 hours and dried by appropriate circulation. Thereafter, the operation of the device is permitted. Nominal data with tolerances are valid after a warm-up time of at least 30 minutes at an ambient temperature of 23°C. A value without tolerances is a guiding value for the device.

Maintenance

The outside of the device should be cleaned regularly with a soft, non-linting cloth. The Diamond Shield should never be cleaned with alcohol or other solvents.

Feature Identification

Front Side

- 1. LED
- 2. LCD Display
- 3. Mute Switch

- 4. START Button
- 5. On-Off Switch
- 6. EXIT Button

- 7. OK Button
- 8. Arrow Keys

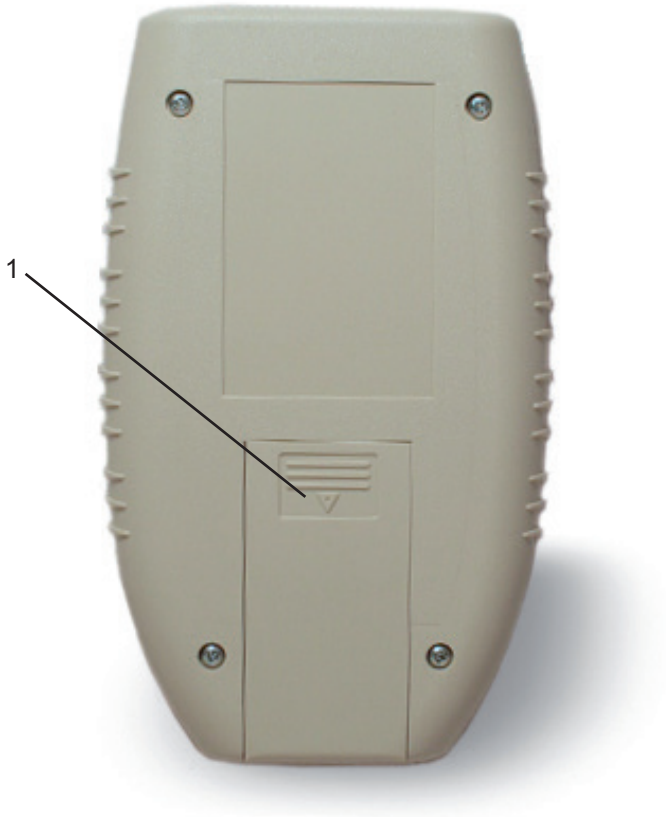
Connection Side

1. Grounding connection
2. Chip Card Slot
3. Connection for electrodes (ie bracelets)

Back Side

1. Battery compartment

Bringing into service

Insert the battery

1. Open the battery compartment cover on the back of the device.
2. Insert the supplied battery. Please keep in mind that the battery has to be inserted correctly. Meaning the plus has to go to plus and minus to minus. (If however the battery is inserted incorrectly, the device will not be damaged by that.)

We recommend using alkaline batteries or rechargeable batteries.

Connecting the cables

1. Connect the connection cable for the bracelet electrodes in the right socket (above LED).
2. Connect the bracelets.
3. If you want to work with grounding: Connect the grounding cable to the left socket. Then connect the ground wire to the outlet.

Operation

Switch On

Pressing briefly the „Power“ button will turn on the device.

Switch Off

Press the „Power“ button and hold it for about 1 second until the device is turned off.

Automatic Switch Off

After 12 Minutes without user or program activity the device switches off automatically to save battery power.

Charge condition of the battery

The state of charge of the battery is displayed on the LED. If the device is switched on, the LED flashes. When the battery is full, the flashing frequency is slow. When the battery is empty the flashing frequency is faster.

If the battery is empty, a message appears, and the device can no longer be used.

Keypad lock

The device has a keypad lock to prevent accidental operation of buttons. To activate or deactivate the keypad lock briefly press the POWER-button and then press OK. The device will display the status of the keypad lock then.

Activate / deactivate the loudspeaker

Press the sound button to switch the loudspeaker on or off. The current state is then displayed briefly.

Menu Navigation

Navigation through the menu structure of the Diamond Shield is simple and identical for all functional areas:

- The currently selected item is marked with a small black rectangle. Use the up arrow key („+“) to select the next higher element and the down („-“) arrow key to select the underlying. If you have already selected the last item of the menu and press the down arrow key, the first entry in the menu is selected. Accordingly, the system jumps to the last entry, if you press the up arrow key when the first entry is already selected. This allows you to switch quickly between the first and the last entry.
- A scroll bar on the right edge display appears when the current list has more entries than can be displayed. The scroll bar then indicates where you are in the list.
- Confirm your selection with the OK button. Then a new display appears that corresponds to the selected item.
- Press the EXIT key to return one level back (if you are not already in the main menu).

Important Terms

- **Frequency program:** A frequency program includes all properties of a single application such as the frequency, intensity or duration.
- **Sequence:** A sequence is a succession of many frequency programs. For example, there is a sequence with the name „Diamond Shield“, which includes the frequency programs „DS1“, „DS2“, ... up to „DS26“.

Menu „Main menu“

A few seconds after switching on the device the main menu will appear. You can always go back to the main menu if you want to leave the submenus by pressing the EXIT button.

The Main menu shows the following entries:

- **Chipcard:** This item is displayed only if a chip card is inserted.
- **Standard:** Under this menu item there are preinstalled sequences and frequency programs available.
- **Configuration:** In this menu various settings can be undertaken.

„Standard“ menu

Under standard you find the pre-installed frequency programs with which the Diamond Shield Zapper is by default supplied with. Most of these programs consist of a sequence of several sub-programs.

The first level of the menu shows the available sequences. When a sequence is selected and then the OK button is pressed, the frequency programs of the selected sequence are shown.

Executing a sequence

Select the desired sequence with the + and - button. Then press the Start button. to run the entire sequence. First the device will test the electrical contact (see „Contact control“). Then all the individual programs of the selected sequence are executed.

Running a single program from a sequence

Select the desired sequence with the + and - button. Press the OK button. The subprograms of this sequence are displayed. Then select the desired sub-program and run it by pressing the start button. Previously the device checks the electrical contact (see „Contact control“).

Contact control

Before a sequence or a single program is executed, the device checks whether the required electrical contact is made. If the contact is not good enough, you will hear a beep every 3 seconds. By pressing the OK button, the frequency program can be executed anyway. However no effect will occur.

If the required contact is made, the unit displays for a short time the resistance in k-Ohms and then automatically begins executing the frequency program.

Display

- In the headline the sequence name and the name of the current single program are shown in turns.
- Below is the remaining execution time displayed (minutes: seconds).
- Below the intensity is displayed.
 - For frequency programs with current control the actual voltage is displayed on the left side. On the right side the selected intensity is displayed in milliamps (1/1000 amps).
 - For frequency programs without current control only the voltage in volts is displayed.

Operating options

- Pressing the Start Button will interrupt the program.
- Repressing the Start Button will play the program from the point of interruption.
- The intensity can be adjusted using the Arrow Keys up and down.
- Pressing the Exit Button will quit the current application.

Menu „Config“

This menu is used to configure the device. The settings made here are retained even after turning off or changing the battery.

- **Language:** Select from English or German. The input is done with the arrow keys left and right.
- **Contrast:** Adjusts the contrast of the LCD display. The input is done with the arrow keys left and right.
- **Const amps:** enables or disables constant current control. If you press the left and right arrow then the input is done. If ‚NO‘ is selected, the intensity of all the frequency programs (including those that are stored on the chip cards) is usually set on the voltage and does not have a constant current control.
- **Speaker:** the tone pitch of the internal speaker. Use the arrow keys left and right to choose 4 or 5 kHz.
- **HW Rev.:** Displays the hardware version of the device.

Playing a chip card

Chip cards are available as accessories containing preprogrammed frequency sequences. Both these diamond chip cards as well as individual chip cards can be played as follows.

1. Insert the chip card in the device.
2. Switch on the device.
3. Select the menu item „chip card“ and press the OK button.
4. A submenu appears, which displays the contents of the chip card. The entire sequence and sub-programs of chip card can be applied similarly as it is described in the chapter „Standard“ menu.

Running the entire sequence of the program: **Start Button**

Selection of single sub-programs from the sequence: **Ok Button**

Technical Data

General

- Operating temperature: +15° C to +30° C
- Storage temperature: -20° C to + 70° C
- Power supply: 9V alkaline (manganese) battery or comparable accumulator.

Connections

- 3.5 mm socket (stereo) for grounding plug
- 3.5 mm socket (stereo) for signal outputs
- Chip card slot

Characteristic of the output signals

- Frequency range 0.5Hz - 2,5 MHz
- Waveform: square wave signal
- Minimum accuracy of the frequencies: 100 ppm
- Output Voltage: 0.1 - 16V (measured from the 0V line)
- Offset: No offset, positive offset, positive offset + 0.3V
- Modulation up to 254 times the output frequency (but 2,5 MHz must not be exceeded)
- Wobble (maximum up to twice the fundamental frequency)
- Voltage control to keep the output current constant (only at frequencies above 20Hz)
- Internal resistance: 1 kOhm

For disposal

The Diamond Shield device must – as well as other electronic devices – not be disposed with normal household waste. You can bring it to the local separate collection system for electrical appliances or return it to the address below (unfree postage will not be accepted)

Please also note that the disposal of batteries is prohibited on the normal household waste for environmental reasons. You can bring them to the public collection places or wherever batteries are sold.

You are welcome to return the battery to the following address (unfree postage will not be accepted):

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© September 2013 DermaVit GmbH & Co. KG
All rights reserved
Any kind of reproduction only with
approval of DermaVit GmbH & Co. KG

DermaVit GmbH & Co. KG
Unterer Anger 15/Rgb.
80331 München
Germany

Tel: +49 (0) 89-18 90 84 6-0
Fax: +49 (0) 89-18 90 84 6-15
E-Mail: info@derma-vit.com
Web: <http://www.derma-vit.com>